


Pipe Surge And Water Hammer Apparatus – C7 MkII

This free-standing accessory to the F1-10 Hydraulics Bench clearly demonstrates the difference between the phenomena of pipe surge and water hammer and how each is created.

MULTIPLE DEMONSTRATIONS & MEASUREMENTS
INCLUDES VIRTUAL OSCILLOSCOPE
DATA LOGGING

F1-10 Hydraulics Bench

The unit is designed for use with the Armfield F1-10 bench

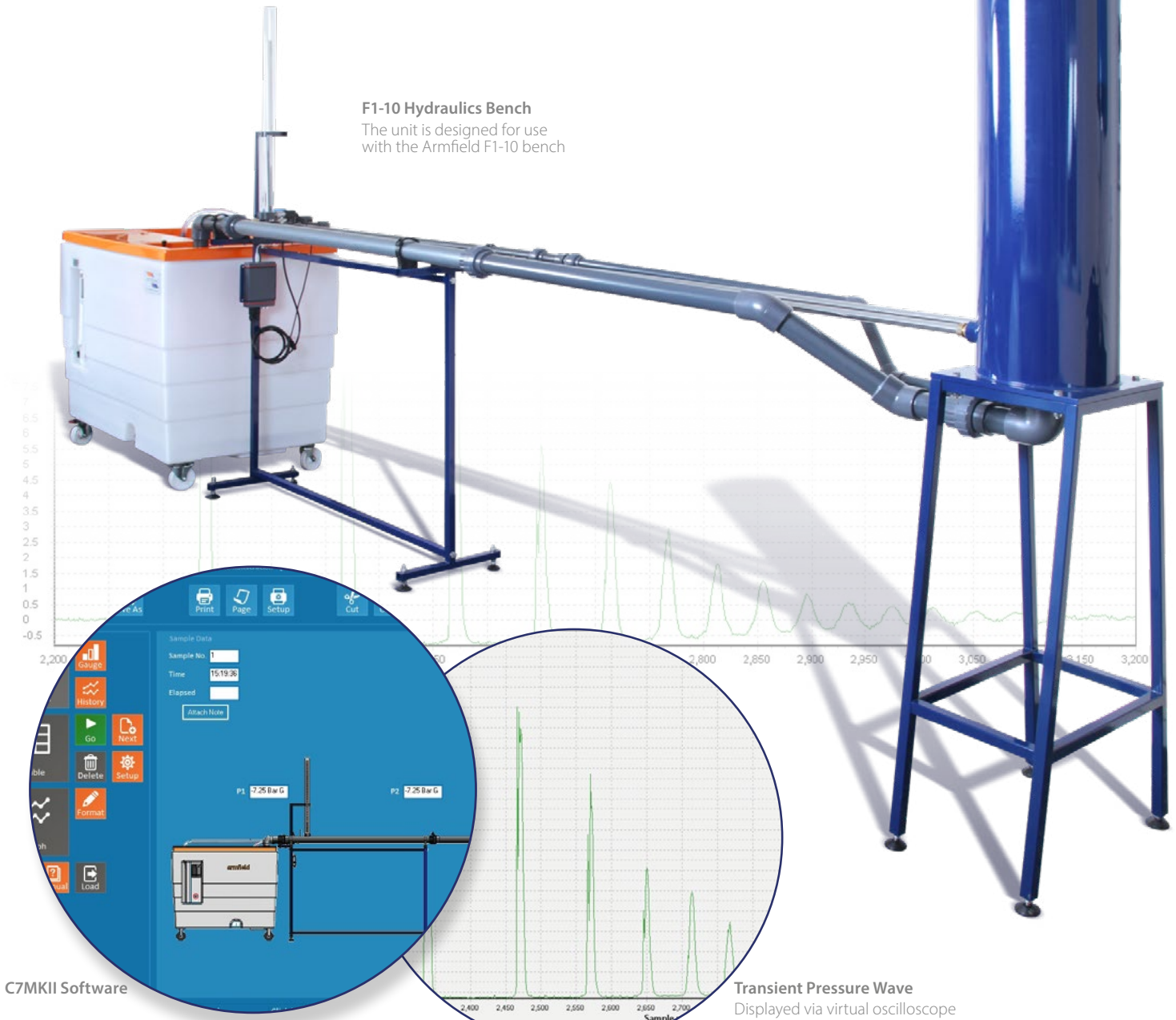

C7MKII Software

Transient Pressure Wave
Displayed via virtual oscilloscope

Experimental content

- ▶ To demonstrate the phenomenon of pipe surge resulting from a change in velocity of the water flowing along a pipe
- ▶ To demonstrate the use of a surge shaft to attenuate any changes in pressure associated with pipe surge and the oscillatory characteristics of water level in a surge shaft
- ▶ To determine the characteristics of water hammer
- ▶ To demonstrate head loss between the reservoir and the surge shaft due to friction in the pipe
- ▶ Interpretation of the traces obtained on the oscilloscope
- ▶ Comparison between theoretical and measured pressure profiles associated with water hammer
- ▶ Determination of the velocity of sound through a fluid in an elastic pipe

UK office - email: sales@armfield.co.uk tel: +44 (0) 1425 478781 (for ROW)

USA office - email: info@armfield.inc tel: +1 (609) 208-2800 (USA only)

Issue: 4

URL: <http://www.armfield.co.uk/c7>

Applications

ChE ME CE IP

We reserve the right to amend these specifications without prior notice. E&OE © 2020 Armfield Ltd. All Rights Reserved

Description

The equipment is free-standing and comprises two stainless steel test pipes connected to a constant head tank with the necessary connections to an F1-10 Hydraulics Bench.

Pipe surge demonstrations are conducted using the first test pipe, which incorporates a transparent surge shaft with scale and a lever operated valve at the discharge end. An additional valve downstream enables the flow through the test pipe to be varied before closing the lever operated valve. An electronic pressure sensor located at the base of the surge shaft enables the low-speed transients in water level to be measured.

Water hammer demonstrations are conducted using the second test pipe, which incorporates a fast-acting valve at the discharge end. An additional valve downstream enables the flow through the test pipe to be varied before closing the fast-acting valve.

The unique fast-acting valve, specifically designed by Armfield, enables water hammer to be generated in a relatively short length of straight

pipe because of the extremely short closure time achieved using a trigger actuator.

Tappings incorporating electronic pressure sensors are located in the test pipe adjacent to the fast-acting valve and halfway along the test pipe. These sensors measure the high-speed pressure transients inside the pipe as the water hammer travels backward and forward along the test pipe. The time delay between the sensors can be used to determine the speed of sound through the water that is attenuated by the elasticity of the metal pipe wall.

All three pressure sensors are connected to a conditioning unit with USB connection for direct connection to a PC (not supplied by Armfield).

The transient pressure waves can be analysed on the PC (not supplied) using the software supplied. The pressure sensors are powered from the USB port on the PC so no additional power supply is necessary.

Requirements


Scale

F1-10 PC USB



- ▶ The user must have a PC with a USB port, running Windows 7 or above.
- ▶ Hydraulics Bench F1-10

Observed oscillations in surge shaft

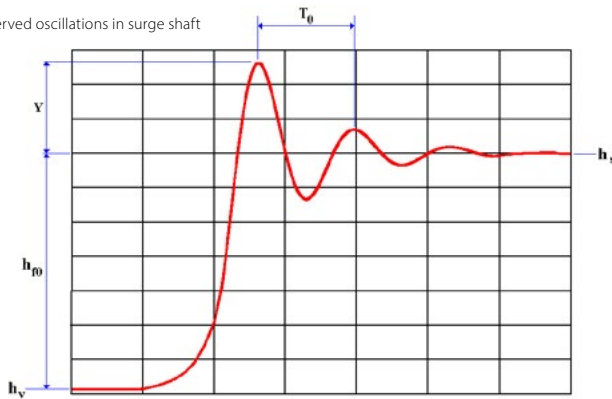

Technical specifications

Max flow rate:	1.6 l/s typical
Max head:	9.0m
Max pump speed:	1,800rpm
Motor power rating:	250W
Pressure sensor:	0 to 15psi
Flow sensor:	0 to 27l/s

Overall dimensions

Dimensions		
Length	3.575m	
Width	0.725m	
Height	1.865m	
Packed and crated shipping specifications		Gross Volume
Gross weight C7MKII	420Kg	5.0m ³
Gross weight C7MKII-10	260Kg	3.5m ³

Ordering specification

- ▶ A free-standing unit designed as an accessory to the Hydraulics Bench to demonstrate the phenomena of pipe surge and water hammer when connected to a Hydraulics Bench
- ▶ Includes two separate stainless steel test pipes, both 3m long, constant head tank, slow-acting valve, fast-acting valve etc
- ▶ A transparent surge shaft (40mm diameter and 800mm high) with scale enables transient water levels to be observed and timed
- ▶ Electronic sensors used to measure pressure transients at two locations in the water hammer test pipe, one adjacent to fast-acting valve and one halfway along the test pipe
- ▶ Data logging software included (requires PC, not supplied)
- ▶ Pressure transients monitored using a PC via a USB connection from the pressure transducers (requires no external electrical supply)
- ▶ Straight metal pipes used, rather than a coiled arrangement, to minimise distortion to the pressure profile

Complimentary products

- ▶ C1MkII: Compressible Flow Bench
- ▶ C3MkII: Multi-pump Test Rig
- ▶ C4MkII: Multi-purpose Flume
- ▶ C6MkII: Basic Fluid Friction Measurements
- ▶ C9: Flow Meter Demonstration Unit
- ▶ C10: Laminar Flow Table
- ▶ C11MkII: Flow in Pipe Networks
- ▶ C15: Computer Controlled Wind Tunnel
- ▶ F1-18: Energy Losses in Pipes
- ▶ F1-21: Flow Meter Demonstration
- ▶ F1-22: Energy Losses in Bends

Ordering codes

- ▶ C7-MkII
- ▶ C7-MkII-10

Armfield standard warranty applies with this product

Knowledge base

- > 28 years expertise in research & development technology
- > 50 years providing engaging engineering teaching equipment

Benefit from our experience, just call or email to discuss your laboratory needs, latest project or application.

An ISO 9001:2015 Company



armfield.co.uk

Aftercare

Installation
Commissioning
Training
Service and maintenance
Support: armfieldassist.com